glo

MANUAL DE INSTRUCCIONES
INSTRUCTIONS MANUAL
MANUEL DE INSTRUCCIONES

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Safety	24
Components	30
Main components	31
Before starting the treatment	32
How to do the treatment	33
Ending the treatment	36
Cleaning and upkeep	36
Customer service and warranty	37

Lippo-massage device

Safety Instructions

Thank you very much for trusting us buying our product; we hope you enjoy your new body care device, an advanced aesthetic treatment for effectively treating cellulite at home, ergonomically designed for intuitive and easy usage, based on the technology used in clinics and medical-aesthetic centres.

In order to provide you with the best possible experience using your new body care device glo910, you can visit our website where you can find all the product information and mode of use in videos. www.glo910.com.

- Read these instructions completely before using your appliance.
- This manual is an integral part of the product. Keep it in a safe place for future reference.

Safety Instructions

- Please, use a swimming cap when use this device for neck and back massage, this appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- This appliance is for domestic use only; it is not for industrial use. Using it for any purpose other than the one indicated could be dangerous. Do not use this appliance as a substitute for appropriate medical care.
- Ask your doctor if you suffer from any health problem, for example, if you use a pace maker.
- If the power cord is damaged, it should be replaced by the manufacturer or by your after-sales service shop or similar qualified personnel to prevent possible hazards.

- Do not attempt to modify or repair the appliance. In the event of any anomaly in the cord or any other part of the appliance, do not use it and take it to an authorised assistance service.
- Do not use the appliance if you are suffering from any skin rash or inflammation.
- Do not use infrared radiation if you suffer from varicose veins, as it could exacerbate the circulatory problem.
- Do not use the appliance in the abdomen area if you are pregnant or believe that you might be pregnant.

- Warning! Do not leave the device on a single point without moving it for more than 30 seconds. The friction generated might cause scratches on the skin or burns in the extreme case of an application on the same point for too long time.
- Do not use parts or accessories not supplied or recommended by the brand.
- Before plugging the appliance into the mains, check the voltage indicated matches that of your home.
- Always unplug the appliance after use and before removing or installing parts and before doing any maintenance or cleaning. Also, unplug it in the case of a power cut. When plugging in or unplugging the appliance, it must always be turned off.
- Do not pull on the cord or hang the appliance from it. Always unplug the appliance by pulling the connection pin, and not the cord.
- Be sure that the cord does not come into contact with sharp edges or hot surfaces. Do not allow the cord to hang over the edge of a table or worktop to avoid children pulling on it and causing the appliance to fall.
- Remove all elements from the inside or outside of the appliance used for protection during transport or for sales promotions, such as paper or plastic bags, plastic sheets, cardboard and stickers.

- Do not use near bathtubs, sinks, showers or other recipients containing water. If it falls into water, do NOT try to pick it up. Unplug it immediately.
- Never leave the appliance outdoors and exposed to the elements.
- Do not touch the appliance with wet hands when it is plugged in.
- Never leave the appliance on a wet surface or cloth while it is plugged in.

RECYCLING

If you decide not to use the grooming kit again, it is best to put it out of service. To do so, unplug the appliance and cut the power cord. In addition, it is best to disable the parts of the appliance that could be dangerous, especially for children, as they might use it improperly during play.

IMPORTANT: When you want to dispose of the appliance, NEVER throw it in the rubbish bin. Take it to the recycling or waste collection centre closest to your home for processing. You will thus be helping to take care of the environment.


This device complies with Low Voltage Directive 2006/95/EC, with Electromagnetic Compatibility Directive 2004/108/EC and Directive 2011/65/EU concerning restrictions on the use of certain hazardous substances in electrical and electronic equipment.

Lippo-massage device
Use

Components

- 1- Main body of the appliance
- 2- On/off massage switch
- 3- On/Off switch with **SLIM tech®** technology (Mod. glo910+)
- 4- **SLIM tech®** technology (Mod. glo910+)
- 5- Variable massage intensity control
- 6- Selected massage intensity luminous indicator
- 7- Placing shaft of the massage rotors
- 8- Button for extracting and securing the massage heads
- 9- Massage heads
 - a. CELLU-ACTIVE (*)
 - b. CELLU-DETOX (mod. glo910+)
 - c. SOFT-MASSAGE (mod. glo910+)
 - d. BODY BRUSH (mod. glo910+)

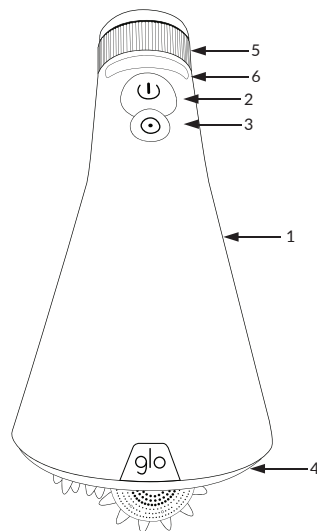
9.a

9.b

9.c

9.d

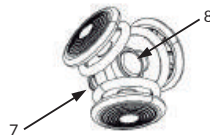Standard (*) glo910 comes with the most versatile and effective treatment body beauty head of the glo 910 series; you can complete and create your own personalized beauty massage treatment acquiring the heads you want separately.

Main components

The **glo 910 series** has been designed to be an **effective therapy to stop the cellulite process**, restoring the correct venous and lymphatic circulation, and eliminating the fat nodules that form the orange peel. Its **powerful engine** performs a rotational movement and strong, but painless, vibration which makes the treatment a pleasant and effective experience.

The glo910 model is equipped with **SLIM tech®** technology. This technology enhances the high efficiency of the fat-reducing massage with a combination of infrared and RED LEDs with different wavelengths (from 630 to 940 nm) that simultaneously penetrate various depths (from 1 to 4 cm). It reaches adipocytes during massage, favours their dissolution and liquefies them to be eliminated through the lymphatic system together with the retained fluids mobilized by massage.

The combination of massage heads, treatment intensity control and thermal energy allows multiple applications:

- **CELLU-ACTIVE:** is ideal for helping to remove cellulite and sagging, stimulate blood flow and return, tone muscles, perform lymphatic drainage, reduce volume, relieve back, sciatica and lumbago pain, muscle relaxation and prewarming, etc. It allows specifically treating abdominal fat to achieve a flat belly effect, thinner waistline and firmer skin.

General Recommendations for Use: Always apply on dry skin. Suitable for legs, buttocks, abdomen and arms.

- **CELLU-DETOX:** Actively stimulates draining and smooths aqueous cellulite. This head is suitable for use with creams and oils, multiplying its results. It can be used for relaxing massages and specific treatments, increasing the penetration rate of these products.

Cellu-Detox is indicated for Athletes as a treatment prior to strenuous exercise, avoiding injuries by pre-warming the muscles to exercise.

This head is also suitable for foot massage (reflexology), which produces a great relaxation and wellness feeling for tired feet with effect on the rest of the body.

General Recommendations for Use: Use with or without cream. Never apply on the back.

- **BODY-BRUSH:** New body scrub head with soft, thin and flexible natural fibre bristles that exert pressure on the body, removing dead cells and impurities to make healthy skin emerge to the surface. The skin regains its youthful, radiant and luminous tone.

Dry brushing tones the body and activates the lymphatic system and blood flow. It has many benefits amongst which fighting against cellulite and dry skin, even at elbows and knees; in addition, it helps preventing varicose veins, stimulates the circulation, helps eliminating toxins and makes the skin look radiant and silky.

General Recommendations for Use: Dry use. Ideally, do this treatment once or twice a week, before bedtime or before the bath/shower. Areas with cellulite: energetic and frequent brushing.

- **SOFT-MASSAGE:** head with soft-foam, flat-surfaced small cylinders, ideal for gentle scrubbing of the skin and cervical massage. In general, for soft and delicate areas, relaxing massages and relief of discomfort in legs with varicose veins. Gently scrubs making the skin light and reaffirmed.

General Recommendations for Use: Use at low speed. Always use dry without creams.

Before starting the treatment

- The skin must be clean and dry
- Choose the suitable head for the desired treatment and place it in the shaft.

PLACEMENT AND REMOVAL OF THE HEAD *(Figure 1)*

- a) To position the head, insert the selected head shaft (7) on the shaft at the base of the device, and secure it by pressing the button (8) located in the centre. Pressing or forcing the anchoring system is needless.
- b) To remove a massage head, press the button (8). Once released, pull it outwards.

Figure 1

How to do the treatment

- Plug in the device into the mains and activate the On-Off button (2) which will light up.
- **Select the massage speed** by turning the knob (5); the selected level is easily seen in the bright viewfinder (6). The head will start rotating. (Figure 2)
- (Model glo910+) If you want to activate the **SLIM tech®** technology (4), press button 4. The button indicating that the thermal power is enabled will light up; in addition, the RED LEDs will be seen through the energy emission area (4) (Figure 3).

RECOMMENDATIONS:

- Apply low-intensity massage in areas where the skin is delicate (neck, breasts, arms...), as well as for the feet, if you have varicose veins or if you are pregnant.
- For the rest of the body, begin with low speed and gradually raise it according to your sensitivity and the head you are using.
- The device is strong enough by itself; excessively pressing it against skin is not necessary; hold it firmly and move it through the skin. You can push harder in some areas such as buttocks.
- At first, bruises may appear with the **CELLU-ACTIVE** (9.a) head, since there may be many entrenched cells; this is normal and gradually improves with treatment. The cause is the intensity applied to deep muscle massage and connective tissue, although it is important to note that not all deep tissue massages end in bruises. If you have capillary fragility and therefore, propensity to present hematomas due to a massage of deep tissue, the smooth Cellu-Detox flat roller (9.b) should be used and begin with a gentle massage intensity.

MESSAGE DURATION:

- The success of this treatment lays on consistency; 15-20 minutes devoted to your care are enough, but you can spend longer time if you prefer. At first, you can massage every day, and later, spread the massages in less weekly sessions, but it depends on you. Look for the ideal moment for a massage along your daily routine; it will help you be constant.
- This device has a built-in security system, so it will automatically shut-off after 30 minutes use. If you wish to continue the treatment, just wait a few minutes (approximately 5) and switch the appliance on again. Continue using normally.

- You can switch the massage heads so that you do not get tired, for example, if you use the **CELLU-ACTIVE** head to fight cellulite, replace it with the **BODY BRUSH** for one day and with **CELLU-DETOX** for another day, increasing the treatment effectiveness.
- The massage should be done in the blood flow direction, in order to favour toxin elimination through the lymphatic system (usually upward). (Figure 4)

Figure 2

Figure 3

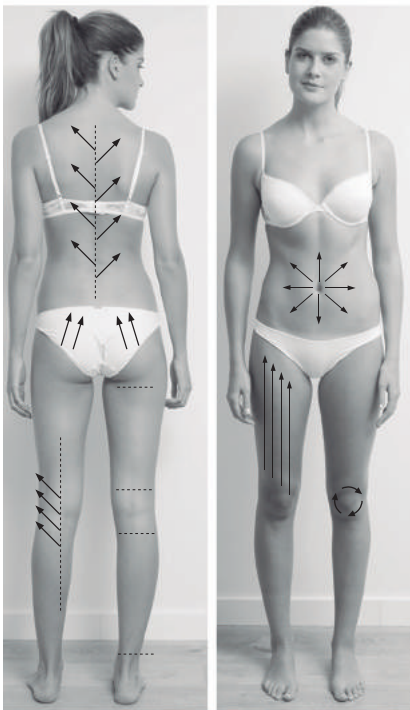
Figure 4


Recommendations for each body area:

! Do not touch bones when massaging; you can hurt yourself.


LEGS: move the device from the ankle towards the knee. From the knee, move the device up to the buttocks and waist.

ARMS: Raise your arm, and move the device from elbow to armpit.

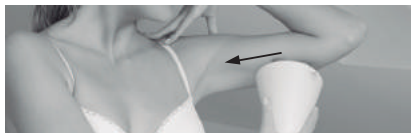
ABDOMEN: in the abdomen, perform the massage at low intensity and slowly move the device from the navel outwards.

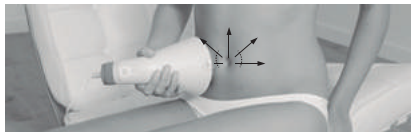APPLICATION PROTOCOL: Discover in www.glo910.com an Aesthetic Dermatology Technique for professional fat reduction results developed in cooperation with a Medical-Beauty Center. This application protocol helps reproducing similar way of massage procedures done at Specialized Centers achieving professional results on the comfort of your home.

BACK: Back massage should be done lying face down; therefore, the assistance of another person is needed. Start massaging in the centre of the lumbar region by moving the device sideways for a few seconds and perform the massage on both sides of the spine up to the neck area.

If you suffer from lumbago, you can use the device to relieve pain. Length, head and strength of this massage will depend on the pain intensity, but it is not advisable to treat this disease for more than 3 min. Then you can apply an anti-inflammatory ointment in the area with the CELLU-DETOX head. If you have the glo910 model, use infrared heads in this treatment.

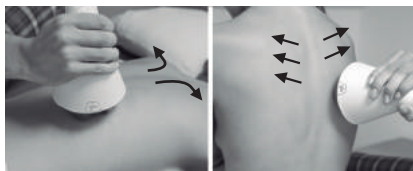
NEVER massage on the spine.

DO NOT leave the device on a single point without moving it for more than 30 seconds. The friction generated might cause scratches on the skin or burns in the extreme case of an application on the same point for too long time.

FOOT MASSAGE: relieves fatigue and stress accumulated in the plant area with a gentle massage at minimum intensity and the CELLU-DETOX head.

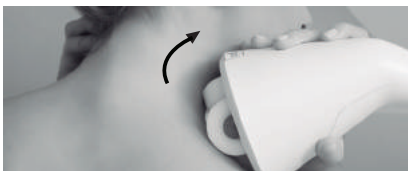Remove hair from the face before applying the device to prevent it from tangling with the head and getting hurt. It is advisable to use a pool or shower cap for more safety and comfort.

To reaffirm the neck, start in the centre of the chin and continue towards the outside of the face, under the jaw, very slowly and without applying pressure.

After intense massage, skin pores are open increasing the effectiveness of treatment creams, so it is the ideal time to apply them. It is even possible to further enhance its effect using the CELLU-DETOX head, for application; its deep massage helps multiplying the effects of your treatment creams and oils. Clean the remains of cream from the head after this treatment.

Ending the treatment

Once the massage is over, turn off the device by pressing the buttons activates at the beginning of the treatment. The massager will stop. Then, unplug it from the mains.

Cleaning and upkeep

- After unplugging the device, use a dry cloth to clean the exterior.
- Do not use solvents, detergents or abrasive products.
- When not in use, keep the device in a cool and dry place.
- To prevent the cable from becoming damaged, do not twist it, pull it or wind it around the appliance. If it is damaged, have it replaced immediately by qualified personnel.

CLEANING THE MASSAGE HEADS:

- Its is recommended to clean the heads with a small brush or a dry cloth. Specifically, they can be washed thoroughly before the next use.

Customer service and warranty :

For further information about glo910 and warranty conditions, please visit our website www.glo910.com.

If your device glo910 + glo910 or is broken, damaged, requires any repair, please contact your local glo910 Customer Service at support@glo910.com.

Guarantee conditions:

- In Member States of the EU, the owner is also granted the rights due to transposition of the 1999/44/EC Directive in the respective national legislation, and, in any other country, by the law in force regarding consumer protection.
- The warranty comes into force on the date indicated on the product invoice or purchase receipt.
- Repair will be carried out completely free of charge, for both materials and labor, during the first two years.
- The holder shall present the product invoice or purchase receipt, unless you register the product at www.glo910.com, so even in the case you lose the invoice, you will be covered by the warranty.
- Damage is not covered by this Guarantee if it is caused by using the equipment in a manner not described in the instruction manual and in general for use other than domestic use, knocks from being dropped, external agents (lime scale, chemical products, water damage etc.).
- Any aesthetic damage caused by use, normal wear and the periodic maintenance operations described in the instruction manual, which may be carried out by the user, is also excluded.
- Supplies (bags, brushes, batteries, bulbs, etc.) are not covered by this warranty.

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